

## One-Step ahead - The obturator without a handle.

One-Step Obturator is the logical next step for the obturator technique without a handle and the use of tweezers making it easier to insert the obturator into the canal.

- A root filling system that is easy to use and gives a perfect seal.
- A product with precise dimensions to make sure it fits the canal.
- A technique with a single insertion that is simpler and more effective than lateral condensation.
- Easy and fast to heat and prepare for a root canal.

**We went even further and designed the One-Step system to make your life easier.**

It is so quick and easy to learn. No longer is there the need to use lateral spreaders and multiple accessory gutta-percha points. No risk of voids.

The ability of One-Step with its excellent flow properties follows the natural canal anatomy and lets you obturate even challenging cases with confidence. When heated, the gutta-percha becomes sticky and tacky, with excellent flow characteristics. It obturates the canal and available lateral and accessory canals easily.

The core is thinner than other leading brands, allowing greater flexibility to navigate curved canals and there is plenty of gutta-percha for denser fills, up to a .10 taper.

Without a handle it is even easier to get an overall view.


The One-Step heater will hold 4 obturators which can be heated at the same cycle. Heating time is 1 minute.

With a universal one-size taper, One-Step obturators have been designed to accommodate most common root canal preparations. It makes it easy to work with and reduces the need for a large stock.


One-Step comes with tip sizes from #20 through #60, matching verifiers, and the One-Step heater.

Unlike other leading carrier brands, One-Step obturators are injection molded to obtain consistent, accurate tolerances required in a root canal treatment and in conformity with ISO 6877.

**With One-Step you are sure to make the perfect apical seal.**

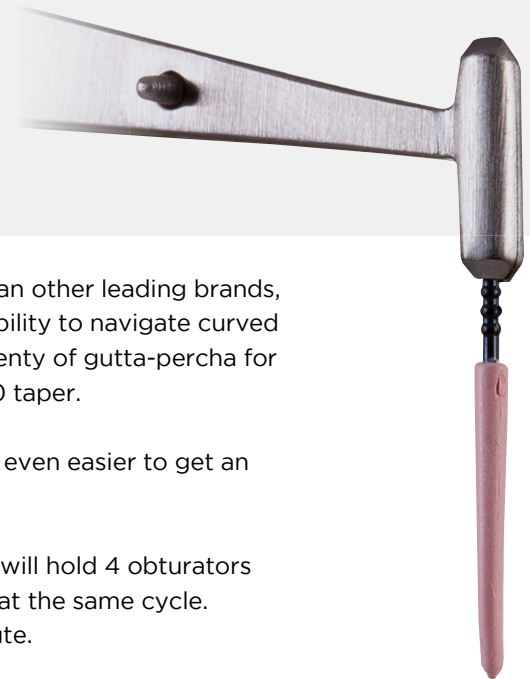

*One-Step obturator*

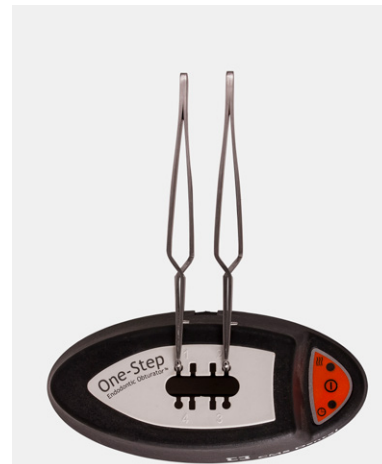

*One-Step heater*



*Root canals prepared for obturation*



*Choosing obturator size*



*Root canals filled with One-Step*



**CMS Dental** has been manufacturing root canal products for satisfied customers worldwide for more than **25 years**.

**Safety, precision, and efficiency are our core values** and should be the focus when dealing with patient care.

The goal has always been to keep it simple, to make treatments safe and dental choices easy. The less technical, the easier it is to work with the products and the less stressful it is for the patients.

The endodontic products are injection molded, which results in precise dimensions and thus makes the dentist's work easier.

The company is using natural gutta-percha and the finest raw materials during manufacturing, which gives the product the properties that ensure a dense and durable root filling.

**CMS Dental is a family-owned company with production facilities in Denmark.**

The company has a top priority on recycling raw materials and has done so since the very beginning 25 years ago. The company, therefore, recycles materials as much as possible without compromising on the quality of its products. CMS Dental aims to achieve greater sustainability in production through further recycling and choosing specialized partners in the local area.



**More information:**

**CMS Dental A/S**

Elmevej 8, Glyngøre

7870 Roslev

Tel. +45 3257 3000

[info@cmsdental.dk](mailto:info@cmsdental.dk)

[cmsdental.com](http://cmsdental.com)